

### Product Information


Stock Number : **AZ0012**  
 Packing Material : Robusta C18 5 $\mu$   
 Length - ID : 250mm x 4.6mm

Serial Number : ~~XXXXXXXX~~  
 Batch Number : ~~XXXXXXXX~~  
 Frit : 2 $\mu$

### Test Conditions

Mobile Phase : Acetonitrile/Water 58/42  
 Flow Rate : 1 mL/min  
 Detector : UV - 254nm  
 QC Station/Tested By : QC-1 / MM

Temperature : Ambient  
 Pressure (psi) : 1932  
 Sensitivity : 1 AUFS  
 Shipping Solvent : Mobile Phase

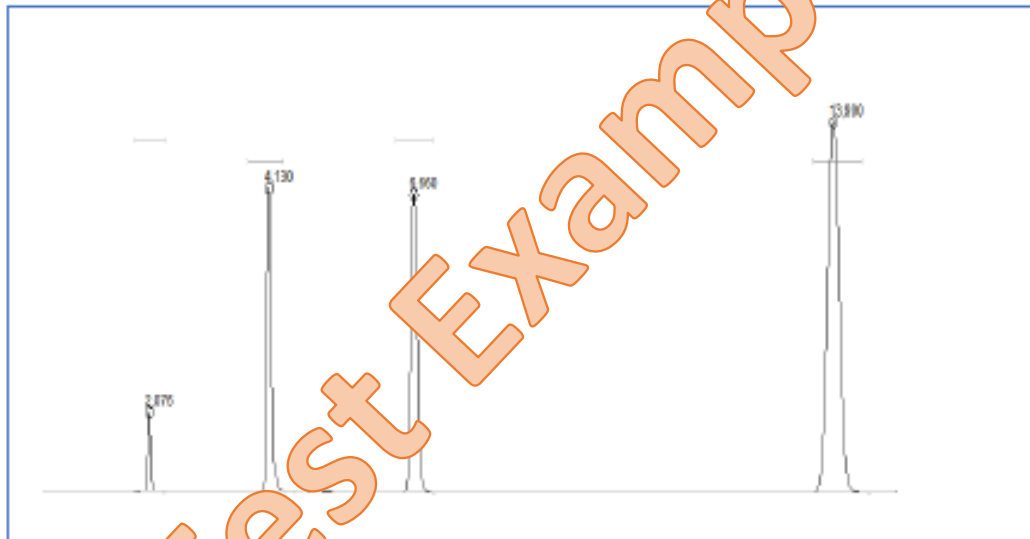

Peak Number	Retention	Capacity Factor	Symmetry	Column Plates	Efficiency (plates/meter)
1	2,08	0,00	1,3	14932	59728
2	4,13	0,99	1,2	19285	77140
3	6,66	2,21	1,2	20308	81232
4	13,99	5,74	1,0	23825	95300

Sample : RP-Mix-D4      Inj. Vol. : 5.0  $\mu$ L

Components	Conc (mg/mL)
1 - Uracil	0.015
2 - Phenol	0.7
3 - N,N-Diethyl-m-Toluamide	0.8
4 - Toluene	4.0